

NEED A MASK SEWING DIRECTIONS

10 Easy Steps

For more information and a video demo go to: www.needamask.life

STEP 1 Cut three sets of the template (two light cotton, two dark cotton, and two polypropylene) and cut 24" of elastic.



STEP 2 Sew each pair together (right sides together for cotton fabric).




STEP 3 Turn right sides of cotton fabric out and top stitch along the edge. Clip the Polypropylene close to the edge.



STEP 4 Open Cotton mask pieces and match front to front. Add the polypropylene to either back side. Sew across top and bottom of mask.




STEP 5 Open a paper clip in the middle and place in the nose section of the polypropylene. Zig zag across to hold in place.

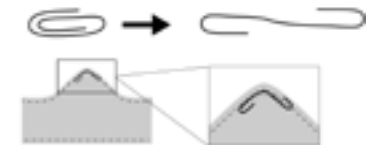

STEP 6 Turn right sides out (keeping polypropylene fabric inside) and top stitch top and bottom of mask. (Leave sides open)



STEP 7 Roll the sides inward 1/4" and top stitch all the way around the mask.

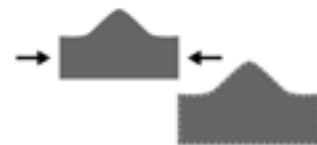

STEP 8 Fold the sides of the mask in 1" and top stitch down to create a casing for elastic.



STEP 9 Feed elastic through casing in a circle with ends hanging on the bottom of the mask.



STEP 10 Tie the bottom of the elastic in an adjustable knot.

